



**CONTINENTAL METAL PRODUCTS**  
Stainless Steel Healthcare Equipment  
Meeting the Demands of the OR

# "QUICK-TEMP"

## Temperature Monitoring System For Blanket and Solution Warming Cabinets

healthcare division

Call us NOW: 781-935-4400

Proudly Made in the  
USA for over 65 years



**NEW**


### "QUICK-TEMP"

**Automatic Temperature Monitoring System Option**  
Blanket and Solution Warmer Temperature Log  
Temperature Log Reporting Made Simple – No More Handwritten Logs

### "QUICK-TEMP"

temperature monitoring system  
now available as an option.  
Requires no additional software.

Download automatically begins  
when a USB flash drive is  
inserted into the front control  
panel. Data is then ready to be  
analyzed on any computer.  
Simple and user friendly.

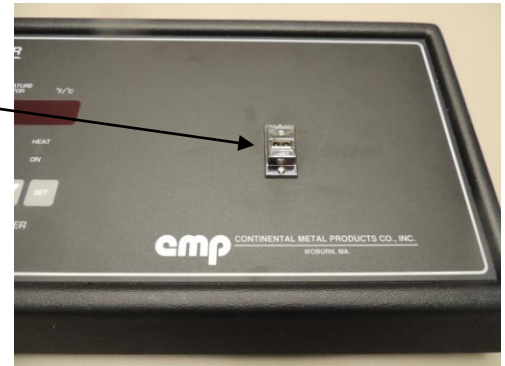

### "QUICK-TEMP" User Friendly Operation

1. Insert USB flash drive into USB port on the front control panel
  2. Temperature log data will automatically download when "LOG" appears in display
  3. Remove USB flash drive
  4. Insert USB drive into any computer to view temperature data
- Alarm conditions are easy to observe and are marked by a large "X" in the report
  - Report contains Corrective Action Taken/Resolution column for clinical staff to make appropriate notes

Action Taken Resolution – provides a space for Clinical personnel to make note that Facilities Management/Plant Services were called and corrective action was taken to meet regulatory guidelines

- Option available on all CMP Blanket and Solution Warmers
- Requires no additional software
- Downloads automatically when standard USB flash drive is inserted into front panel
- Data is ready to be analyzed on ANY computer
- ½ hour sample rate
- 1 year temperature log history
- CSV file for easy viewing on standard spreadsheet software
- All warmer temperature data and set points are displayed in report
- Report displays if an alarm condition occurs
- Area to note action steps taken in alarm situation
- Assists in meeting Blanket Warmer Guidelines for Temperature Monitoring and Regulatory Requirements\*

### Sample Report

Continental Metal Products – Blanket/Solution Warmer Temperature Log										
Date	Time	Upper Comp. Temp.	Lower Comp. Temp.	°C/°F	Upper Comp. Setpt.	Lower Comp. Setpt.	Upper Comp. Limit	Lower Comp. Limit	Alarm Condition (X)	Corrective Action Taken Resolution
10/12/2012	08:32	105	125	F	105	125	180	180		
10/12/2012	09:02	105	125	F	105	125	180	180		
10/12/2012	09:32	105	125	F	105	125	180	180		
10/12/2012	10:02	105	125	F	105	125	180	180		
10/12/2012	10:32	105	125	F	105	125	180	180		
10/12/2012	11:02	105	125	F	105	125	180	180		

**Continental Metal Products Co., Inc. – 35 Olympia Avenue, Woburn, MA 01801**  
Phone (781) 935-4400 FAX (781) 935-4404 [www.continentalmetal.com](http://www.continentalmetal.com) [sales@continentalmetal.com](mailto:sales@continentalmetal.com)

\* "ECRI Institute recommends that temperature settings on blanket warming cabinets be limited to 130°F (54°C) and that solution warming cabinets be limited to 110°F (43°C). Temperatures above this level unnecessarily increase the risk of burns while providing no added clinical benefit."



**CONTINENTAL METAL PRODUCTS**  
Stainless Steel Healthcare Equipment  
Meeting the Demands of the OR

# Blanket and Solution Warming Cabinets

healthcare division

Call us NOW: 781-935-4400

Proudly Made in the  
USA for over 65 years



**NEW**

## "QUICK-TEMP"

**Automatic Temperature Monitoring System Option**

Blanket and Solution Warmer Temperature Log

Temperature Log Reporting Made Simple – No More Handwritten Logs

## "QUICK-TEMP"

Option available on all CMP  
Blanket and Solution Warmers



CMP has been building high quality Stainless Steel Warming Cabinets for over 65 years and has become the industry standard for hospital and OR blanket and solution storage.

**CMP is happy to configure your warming cabinets to meet your unique requirements**

**CMP is so confident in our quality and craftsmanship that we offer a two year labor and six year parts warranty**

- Stainless Steel Construction
- Electric Heated Compartments
- Over Temperature Alarm
- Separately Controlled Chambers
- Temperature Adjusting Lockout
- Optional "Quick-Temp" USB temperature monitoring system
- Many Standard Sizes Available
- Microprocessor Control System
- Customizable Temperature Settings/Limits
- Wire Shelves to maximize airflow
- Removable Pilaster Strips for easy cleaning
- Optional Wi-Fi and Ethernet connectivity
- Optional 3/4" thick Thermopane Glass Doors
- Custom Sizes
- Digital Display
- UL/CUL Listed
- ± 2°F Accuracy
- Uniform Heating
- Optional Casters
- Made in the USA

**Continental Metal Products Co., Inc. – 35 Olympia Avenue, Woburn, MA 01801**  
Phone (781) 935-4400 FAX (781) 935-4404 [www.continentalmetal.com](http://www.continentalmetal.com) [sales@continentalmetal.com](mailto:sales@continentalmetal.com)